



Save • Someone's **life**



AED-Defi.com
Medical Equipment Manufacturer

Defi[®] 5ST
AED Trainer



Start the device

1

Turn on the power.
The voice instruction starts immediately for guide on how to operate AED.



Attach the pads

2

After you connect the pads,
the trainer starts working with the
selected mode (one of four).

Selection mode

3

Mode 1: Shock indication
and CPR voice instruction

Mode 2: No shock indication
and CPR voice instruction

Mode 3: Shock indication
and no CPR voice instruction

Mode 4: No shock indication
and no CPR voice instruction



Press the button

4

When the trainer asks you to give a shock
(the trainer does not have an electric shock).
Press the red shock button



Advantage

15 Languages selectable	Child Mode	Adult Mode	Program 4 training modes	Sound volume 3 Levels	Carry Case for AED and first Aid

- 15 languages to choose from (by end user)
- 3 volume levels (by the end user)
- Button for adult/child mode
- 4 training modes
- Use of the trainer in the pocket (does not need to be opened)
- Extra compartment in the bag for first aid products
- Extensive voice prompts for the operator
- Manually switchable language
- Runs on standard AA x 4 batteries


■ More than 45 languages are available

**Speaks
your
language**



Technical Specifications

Mode	Adult, Child
Dimensions	Depth (D): 80mm Width (W): 240mm Length (L): 300mm
Weight (with battery)	1.9Kg (4.2 Pounds)

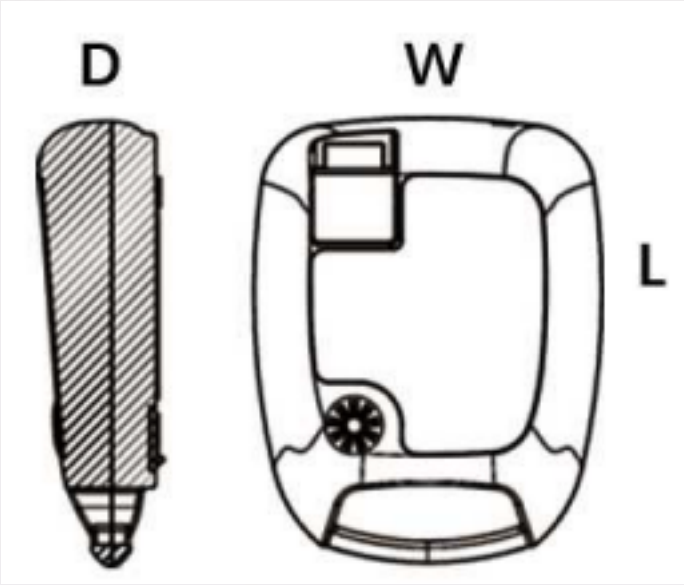


Accessories

■ Battery

4 x AA Battery (6Volt)

■ Reusable pads

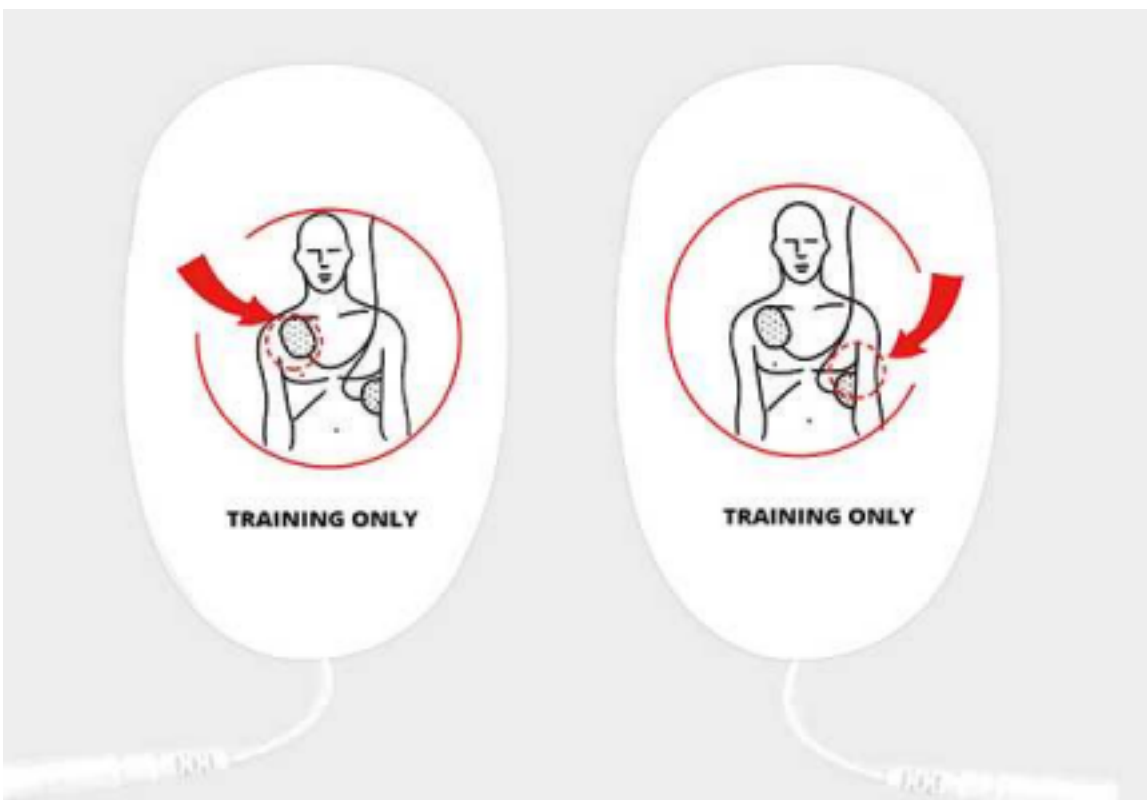

■ Carry Case

- 1 Reflective sign
- 2 Soft Hand
- 3 Speaker protection
- 4 Transparent panel
- 5 Socket Protector

