



Hyperbaric Oxygen Therapy

Henry's Law:

“The amount of dissolved gas in a liquid is proportional to its partial pressure above the liquid.”

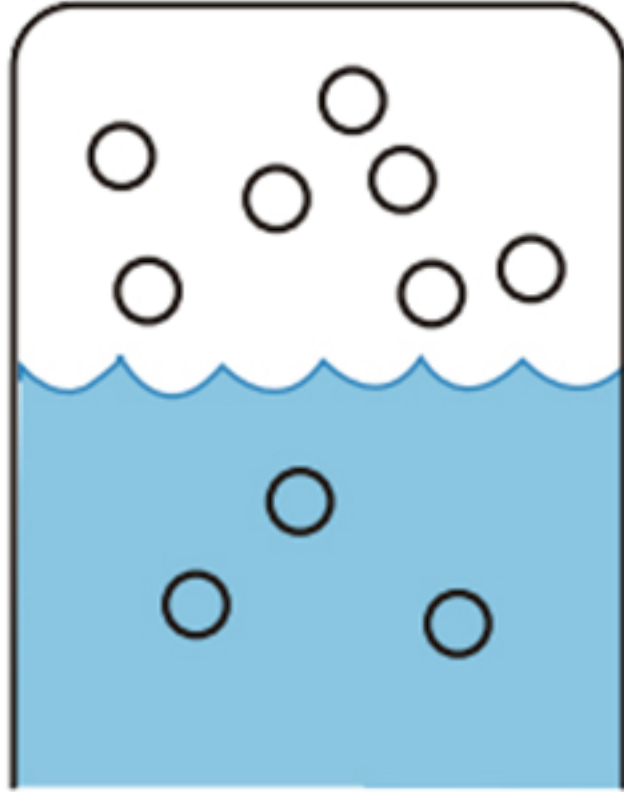

The state of gas dissolved in liquid at one ATA pressure.

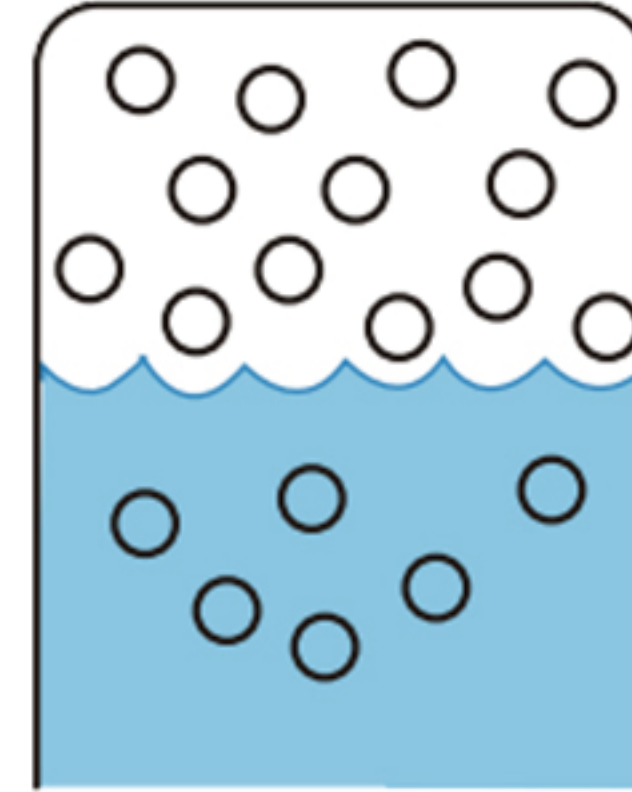

The state of gas dissolved in liquid at two ATA pressure.

Hyperbaric Oxygen Therapy is the medical use of oxygen at an increased level of pressure, (hyper) means increased, (baric) means pressure. Hyperbaric Oxygen Therapy is a treatment which enhances the body's natural healing process by providing an environment which allows the body to absorb much higher amounts of oxygen than possible at normal atmospheric pressure. It increases the oxygen transport capacity in the blood. This noninvasive therapy is the most trusted way to increase oxygen levels to all organs of the body.

Applications of hyperbaric chambers



Adjuvant treatment of some diseases



Rapid recovery after exercise




Family health management



Beauty salon anti-aging

Product Information

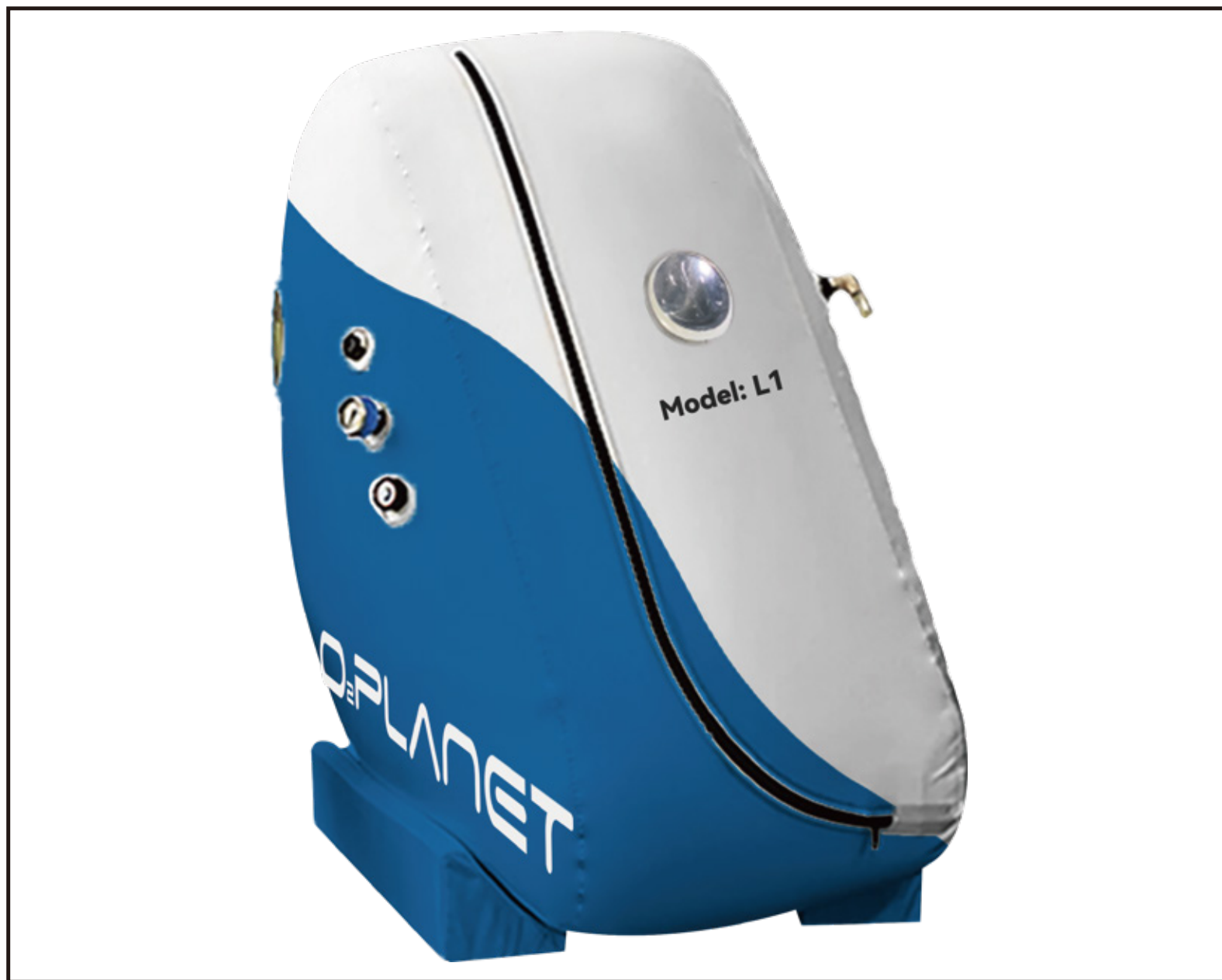

Medium size lying type for 1 person

Chamber Parameters

- Size(L*W*H): 140*100*160cm / 56*40*63"
- Application: Clinic, Gym, Spa, Homecare
- Pressure: 1.3/1.4/1.5ATA (30/40/50kPa)
- Package dimension: 87.5*57.5*47cm
- Net weight: 28kg
- Gross weight: 30kg



Oxygen Concentrator

Imported grade high performance molecular sieve Effectively filter 99% of bacteria Provide you with ultra-pure forest-like oxygen.

- Size(L*W*H): 35*10*65cm / 14*4*26"
- Flow: 1-5L/min (adjustable)
- Power: 500W



Air Cooler

Sufficient power, built-in filter cotton to purify the compressed air, equipped with a silencer to minimise noise.

- Size(L*W*H): 18*12*35cm / 7.2*5*14"
- Power: 200W



Air Compressor


Sufficient power, built-in filter cotton to purify the compressed air, equipped with a silencer to minimise noise.

- Size(L*W*H): 39*24*26cm / 16*10*10"
- Flow: 72L/Min
- Power: 550W

L1 Features

- Elongated L-shaped zippers for easier access.
- Comfortable and safe ergonomic design.
- Multiple vision windows to check what's going on inside/outside.
- Oxygen headset/mask for breathing in high-purity oxygen.
- Two automatic pressure relief valves to adjust the pressure inside and internal/external pressure gauges to monitor the pressure inside.
- Emergency pressure relief valves for fast exit; pressure relief valve operable both inside/outside.
- An air compressor, an oxygen concentrator and an air cooler for the control system.
- One piece molding, user-friendly design, easy to install and operate.

Product Superiority

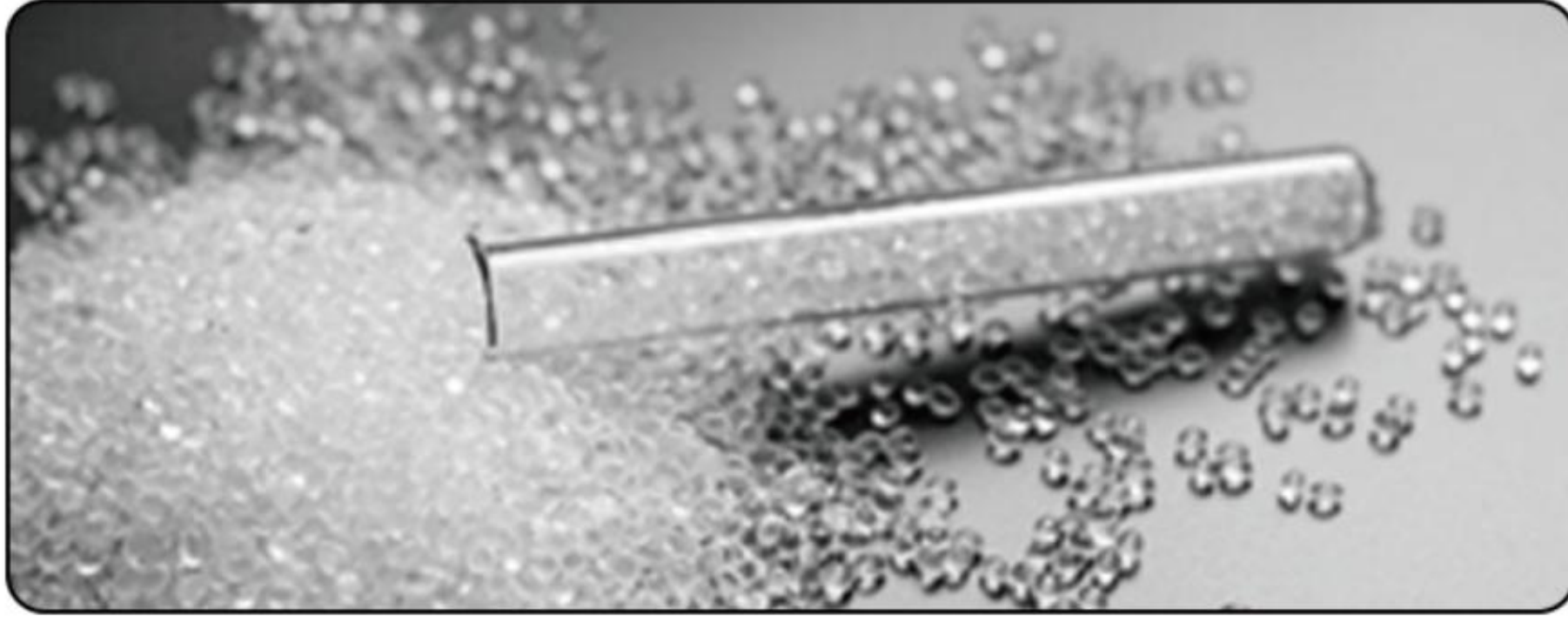

Cabin raw materials

TPU material passed SGS and Rohs environmental certification, high strength non-toxic environmentally Good elasticity and abrasion resistance, tear resistance, high tensile strength.

Welding

Energy-saving and environmentally friendly high-frequency welding a molding process safe and more stable.



Optional outer cover

Half wrap and full wrap cloth covers are available for choose Different Color can be chosen 100% cotton stretch fabric. High wear-resistant and dust-resistant polyamide fiber cloth cover

Zipper

Yoshida Kogyo Kabushikigaisha The world's top YKK brand smooth zipper The originator of the zipper industry, representing the industry standard Adopt Japanese precise process.



3 types of safety valves ensure safe and effective use



Pressure Gauges



Automatic Pressurize Relief Valves



Emergency Relief Valve



Air Deflate Value

[CONTACT US](#)

Product Pictures & Customer Feedback

